

Cleaning Instructions

XLR8 / XLR8+



// CAUTIONARY NOTE

- To help reduce risk of infection and contact with blood and body fluids, use personal protective equipment (PPE) such as medical procedure gloves
- When disinfecting the units, extreme caution must be taken to ensure that no excess fluid is allowed to pool in the areas surrounding the screen or the outside edges of the rubber seal.
- Excessive pressure should not be applied to the screen.
- Be careful not to allow any moisture to enter the pump while cleaning.
- Do not use abrasive wipes on the screen.

// CLEANING PROCESS

1. If there is a Carrying Bag, please remove the pump from the bag and discard.

2. Use Sporicidin® Disinfectant Towelettes/Wipes or any FDA approved disinfectant cleaning solutions/ wipes.


If unsure, a Genadyne Bio-Med representative can help the facility determine the usability of the specific disinfectant.

If necessary, wring out the wipe before using to reduce liquid in wipe.

3. Clean all organic material (visible soil or body secretions) from the therapy unit prior to disinfection.

4. Using the wipe, start at the canister side of the unit.

Wipe down the sides of the unit and then wipe down the front of the unit. Gently clean smudge marks off the area around the screen. Tilt the unit so fluid will not pool onto screen area.

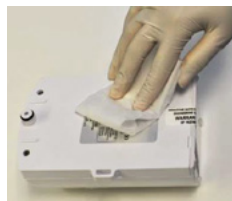

-
5. Wipe down the screen. Extreme caution must be taken when wiping down the screen to not allow fluid into the screen. It should be done very lightly to avoid damaging the screen.
-

6. Lay unit flat with screen area up and allow to air dry.



-
7. After the unit has dried, wipe down the screen with “micro-fiber” cloth to reduce liquid collection on screen and clean off any smudge marks on the screen.
-

8. Look for obvious damage to the housing (cracked/broken), screen, and power port. If there is damage then tag it and follow local procedures for non-conforming products or contact Genadyne for a Return Merchandise Authorization.
-

9. Use wipe to clean adapter, cord, and crevices of the power connector. Wrap the power cord in a looping fashion.



-
10. When the cleaning process is done, dispose the gloves and wipes. Please be sure to wash your hands.
-

11. Move clean unit to another area to avoid cross-contamination. Do not handle clean products with dirty gloves.
-