

# UNO

## Indications

Genadyne UNO is indicated for use in patients who would benefit from negative pressure wound therapy particularly as the device may promote wound healing by the removal of low to moderate exudate and infectious material.

Appropriate wound types include:

- Chronic
- Acute
- Traumatic
- Subacute and dehiscent wounds
- Partial-thickness burns
- Ulcers (such as diabetic or pressure)
- Flaps and grafts
- Closed surgical incision

Genadyne UNO is a single patient use device.

## Cost Effectiveness

The 7 day therapy options provided by the UNO System is a convenient cost effective alternative to other disposable systems. With a longer pump life, an added canister and the option to order disposables separate from the pump, treatment cost is significantly reduced.

Genadyne Biotechnologies  
16 Midland Avenue  
Hicksville  
New York  
11801  
USA

Phone: +1 516 487 8787  
Toll free: +1 800 208 2025  
Fax: +1 877 487 7878  
Email: [info@genadyne.com](mailto:info@genadyne.com)

## Kit Presentation

**Each System comes with everything in one box**

- 1 Pump
- 1 70ml Canister
- 2 Dressings
- 6 Security strips
- Carrying bag
- Power adapter
- USB cable

Everything you need in one convenient box ensures all materials are present when initiating and continuing therapy. Additional dressings and canisters may be ordered separately, guaranteeing that the pump can be used for the full 7 days of therapy.

### Comes in the following range of sizes:

| Description             | Ref       |
|-------------------------|-----------|
| UNO 7 Therapy Kit 10x10 | U-S1010-4 |
| UNO 7 Therapy Kit 10x20 | U-S1020-4 |
| UNO 7 Therapy Kit 10x30 | U-S1030-4 |
| UNO 7 Therapy Kit 10x40 | U-S1040-4 |
| UNO 7 Therapy Kit 15x15 | U-S1515-4 |
| UNO 7 Therapy Kit 15x20 | U-S1520-4 |
| UNO 7 Therapy Kit 15x30 | U-S1530-4 |
| UNO 7 Therapy Kit 20x20 | U-S2020-4 |
| UNO 7 Therapy Kit 20x25 | U-S2025-4 |
| UNO 7 Therapy Kit 25x25 | U-S2525-4 |



Single patient disposable fully functioning  
Negative Pressure System

# UNO single patient disposable fully functioning Negative Pressure System



## The Pump

Lightweight – only 230 grams (with an empty canister).  
 Small – 7.6cm x 11.1cm x 5.7cm – easily and discreetly carried by the patient.  
 Pump cover – for additional discretion and ease of wearability.

7 day Pump options available.

### Therapy Treatment Time

Each pump provides 7 days of effective therapy plus one ‘grace day’.

### 2 Modes of Operation


Continuous – at either 125mmHg or 80mmHg.  
 Variable – at either 80mmHg/30mmHg or 125mmHg/30mmHg.  
 Simple – press just 3 buttons to select and start the required therapy the clinician desires.

### 3 Alert Settings

Leak, canister full/blockage and battery low - you can be confident that the UNO System is providing optimal therapy at all times, with no ‘downtime’ where the system is not operating effectively.

Intermittent pressure levels of other brands = 125-0mmHg (spike)

Variable pressure levels of UNO = 80-30mmHg (gradual)

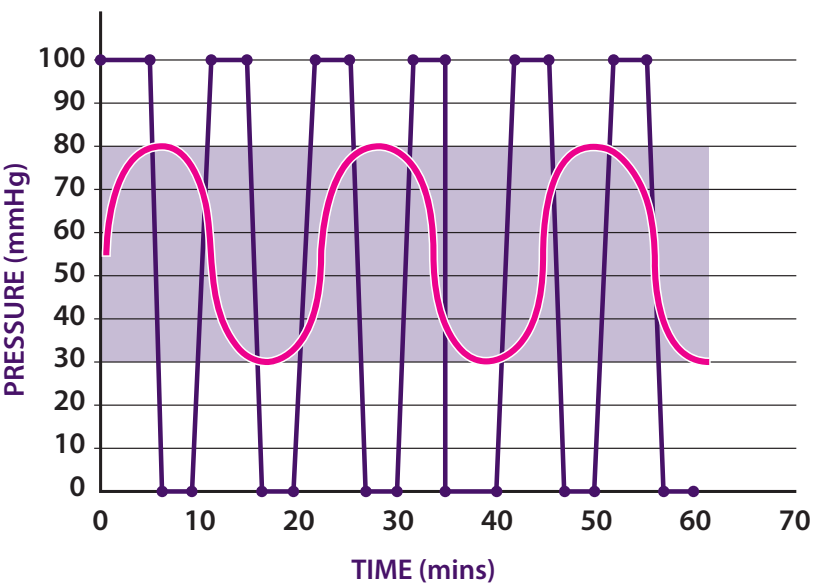

Variable pressure has been proven to produce the best healing rates of any therapy setting

## Canister Presentation

As well as an absorbent dressing regimen, the UNO System also has an integrated 70ml canister which enables the system to easily manage exudate on both low and moderately exuding wounds. It does not rely solely on moisture vapor transmission to manage exudate.

Each canister also includes a gelling and deodorizing agent.



Polyurethane Film

Super absorbent layer

Foam layer

Silicone gel adhesive layer

## The Port Dressing

Each kit comes with 2 individually packed and sterilized dressings, enough for 7 days of therapy, depending upon which option you choose. It is recommended to change the dressing every 3 days or up to 7 days as per physician recommendation. Comes pre-assembled for quick and easy application. Just one quick luer-lock connection to the pump.

**Foam layer** for quick absorption and maintaining a moist environment at the wound site.

**Super absorbent layer** for additional absorption and aids the movement of exudate into the port.

**Silicone wound contact layer**  
 Soft and gentle application and removal for enhanced patient comfort. Can be repositioned to guarantee a good seal on the dressing (no ridges/tracks where exudate could pass through).

**Film layer** allows small amounts of air to pass through, needed for exudate to be pulled through into the canister. Transpires exudate away from the wound.

**Security strips** – each dressing comes with 6 security strips for additional security (if required) in those areas where dressing adherence is difficult for additional peace of mind.

## Port Position

UNO dressings, unlike non-canister systems, can be positioned either above or below the wound site, and will still work effectively.

Non-canister systems need to be changed when exudate approaches the port position.

With UNO, exudate movement into the port is encouraged.

Because of the inclusion of a canister, the dressings DO NOT NEED TO BE CHANGED WHEN EXUDATE REACHES THE PORT.

## UNO vs Competitor: Cost effectiveness and dressing changes

