

HOW IT WORKS

The Uno System provides lateral tension to the wound site or incision line. This force provides downward pressure which may aid in decreasing incidence of edema and hematoma formation and help reduce risk of wound dehiscence.

The Uno System creates a negative pressure environment that improves blood circulation to the wound and incision site. This enhanced circulation facilitates the delivery of essential cells needed for the formation of cellular collagen, elastin, and the extracellular matrix. These components are crucial for closing incisions and promoting wound healing.

INDICATIONS

Suitable wound types encompass:

- Chronic
- Acute
- Traumatic
- Subacute
- Dehisced Wounds
- Partial-thickness burns
- Ulcers (indicating diabetic or pressure ulcers)
- Flaps and grafts
- Closed surgical incisions



BILLING

Home Health

HCPCS CODE: A9272

Beginning January 2024, Medicare will make a separate payment for HCPC A9272 for a disposable Negative Pressure device along with Home Health Prospective Payment 3*.

*Please note that Medicare policies and coverage criteria are subject to change. It is advisable to consult with a healthcare professional or Medicare representative for the latest information on Medicare coverage for wound therapy devices.

Physician (Office) Hospital, Outpatient Dept (HOPD, i.e., WCC)

CPT CODES: 97607

wounds less than 50 sq cm

CPT CODES: 97608

wounds greater than 50 sq cm

Negative pressure wound therapy, (e.g., vacuum assisted drainage collection), utilizing disposable, non-durable medical equipment including provision of exudate management collection system, topical application(s), wound assessment, and instructions for ongoing care, per session;



Genadyne Biotechnologies, Inc.

16 Midland Ave., Hicksville, NY 11801
Tel: +1.516.487.8787 | +1.800.208.2025

www.genadyne.com

If you have questions about our product or need support, or interested in becoming a distributor, please email or call:

info@genadyne.com | +1.516.487.8787 | +1.800.208.2025




A single use, cost effective, 30 day disposable Negative Pressure Wound Therapy System designed for low to moderate severity wounds, and incision management.

the power to heal

THE PUMP

- **Lightweight**
.5 lbs or .23kg
- **Rechargeable**
USB Charger
- **Compact**
7.6cm x 11.1cm x 5.7cm
- **Easy and Quick to Program**
2 Modes of Operation available
Continuous 125MMHg or 80mmHG
Variable 80mmHG/30mmHG or 125mmHg/30mmHg
- **Discreet**
The device can be conveniently and discretely carried by the patient
A pump case is provided for added discretion and enhanced ease of wearability
- **Extended Therapy Time**
The pump provides 30 days of effective therapy

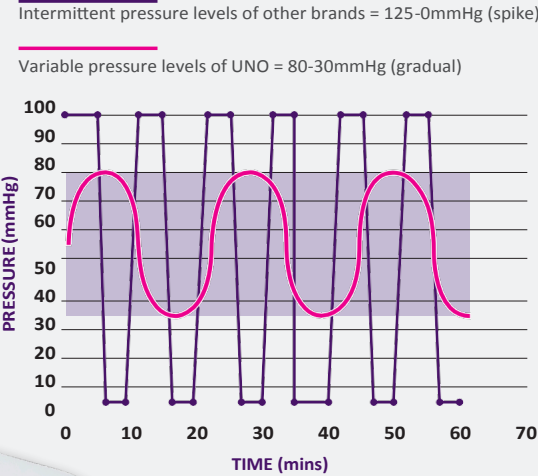

- **Alert Settings**
- Leak
 - Blockage
 - Canister Full
 - Battery Low



THE CANISTER

- **Replaceable Canister**
70ml replaceable canister
- **Exudate Control**
Unlike relying solely on moisture vapor transmission, the system employs this canister for efficient exudate control of low to moderate exuding wounds
- **Enhanced Performance**
Equipped with gelling agent



THE PORT

- **Port Position**
Unlike non canister systems that require changing when the exudate approaches the port, the Uno, with its 70 ml canister, encourages exudate movement into the port and then the canister, allowing for more efficient and continuous operation
- **No Cutting**
Comes pre-assembled for swift and straightforward application
- **Single Luer-lock**
Connection to the canister tube allows the patient to initiate therapy promptly

THE UNO VS. COMPETITOR

	A	B	C
COMPETITOR	GOOD	CHANGE	CHANGE
UNO	GOOD	GOOD	CHANGE

THE DRESSING

- **Apply with ease** Simply peel off the bottom silicone layer and place over the wound
- **Skin Friendly** Silicone layer minimizes damage to skin damage
- **Repositionable** Can be repositioned if desired Comprised of 4 Layers



THE KITS

- **UNO**
- **ADHESIVE STRIPS**
- **DRESSINGS**
- **CANISTERS**

A range of sizes of dressing to handle most postoperative wounds.



PRODUCT CODE	PRODUCT DESCRIPTION	KIT TRAYS	FOAM KIT	70ML CANISTER	CARRY CASE	USB CHARGER	PUMP
U-S1010-6	UNO 30 THERAPY KIT 10X10	5 PER SET	1 SET	2	1	1	1
U-S1020-6	UNO 30 THERAPY KIT 10X20	5 PER SET	1 SET	2	1	1	1
U-S1030-6	UNO 30 THERAPY KIT 10X30	5 PER SET	1 SET	2	1	1	1
U-S1040-6	UNO 30 THERAPY KIT 10X40	5 PER SET	1 SET	2	1	1	1
U-S1515-6	UNO 30 THERAPY KIT 15X15	5 PER SET	1 SET	2	1	1	1
U-S1520-6	UNO 30 THERAPY KIT 15X20	5 PER SET	1 SET	2	1	1	1
U-S1530-6	UNO 30 THERAPY KIT 15X30	5 PER SET	1 SET	2	1	1	1
U-S2020-6	UNO 30 THERAPY KIT 20X20	5 PER SET	1 SET	2	1	1	1
U-S2025-6	UNO 30 THERAPY KIT 20X25	5 PER SET	1 SET	2	1	1	1
U-S2025-6	UNO 30 THERAPY KIT 25X25	5 PER SET	1 SET	2	1	1	1
PRODUCT CODE	PRODUCT DESCRIPTION	KIT TRAYS	FOAM KIT	70ML CANISTER	CARRY CASE	USB CHARGER	PUMP
U-S1010-6S	UNO 30 THERAPY KIT 10x10	2 PER SET	1 SET	1	1	1	1
U-S1020-6S	UNO 30 THERAPY KIT 10X20	2 PER SET	1 SET	1	1	1	1
U-S1030-6S	UNO 30 THERAPY KIT 10X30	2 PER SET	1 SET	1	1	1	1
U-S1040-6S	UNO 30 THERAPY KIT 10X40	2 PER SET	1 SET	1	1	1	1
U-S1515-6S	UNO 30 THERAPY KIT 15X15	2 PER SET	1 SET	1	1	1	1
U-S1520-6S	UNO 30 THERAPY KIT 15X20	2 PER SET	1 SET	1	1	1	1
U-S1530-6S	UNO 30 THERAPY KIT 15X30	2 PER SET	1 SET	1	1	1	1
U-S2020-6S	UNO 30 THERAPY KIT 20X20	2 PER SET	1 SET	1	1	1	1
U-S2025-6S	UNO 30 THERAPY KIT 20X25	2 PER SET	1 SET	1	1	1	1
U-S2025-6S	UNO 30 THERAPY KIT 25X25	5 PER SET	1 SET	1	1	1	1

*disposables also sold seperately