


QUICK REFERENCE GUIDE

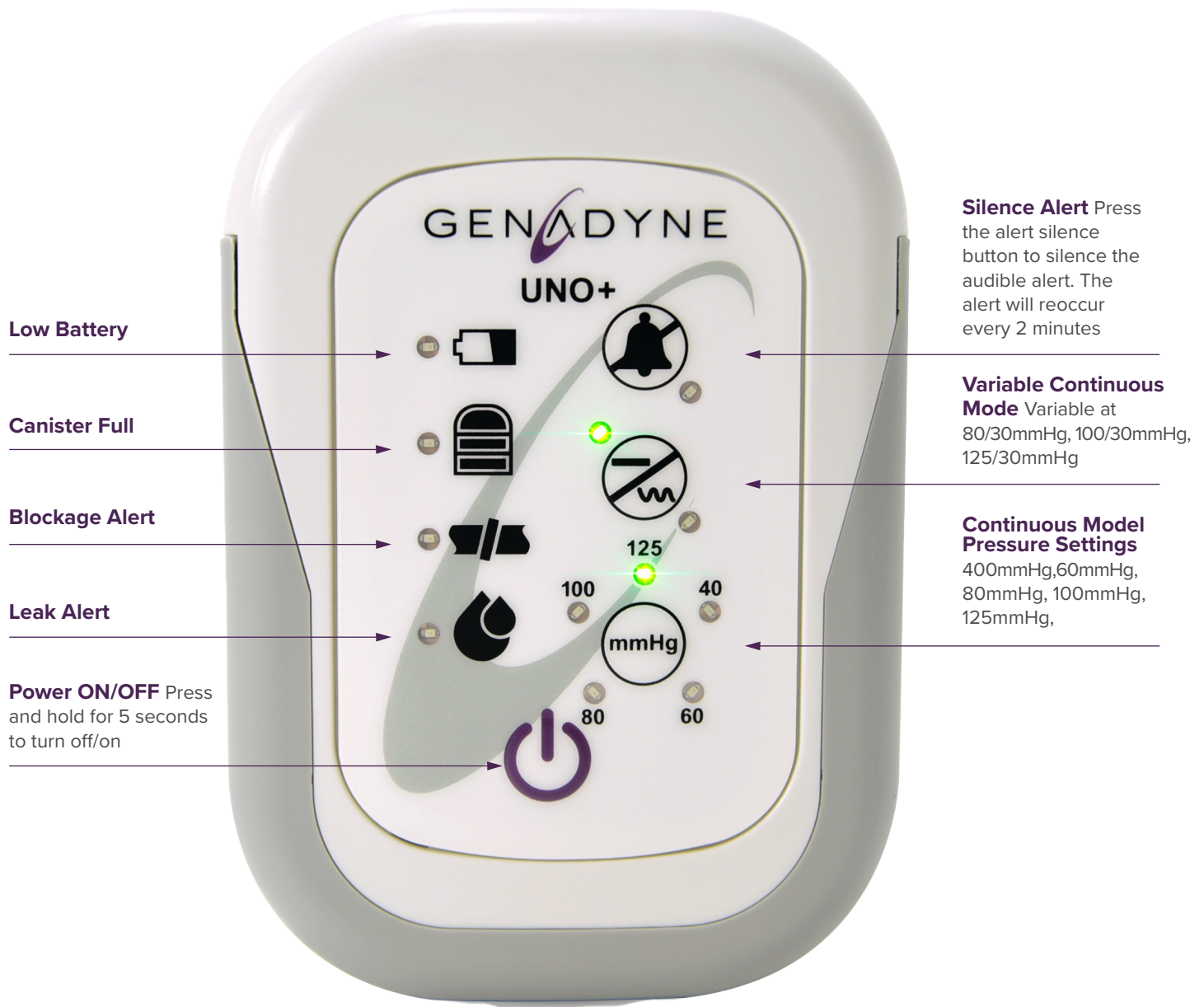

QUESTIONS?

Call: +1.800.208.2025 or +1.516.487.8787

Email: info@genadyne.com | www.genadyne.com

DMR-07-096 rev A

the power to heal

UNO+

QUICK REFERENCE GUIDE

GENADYNE

BUTTONOLOGY



Turning the UNO+ ON/OFF

Press and hold the **Power Button** 3 seconds to power **ON/OFF** the UNO+ pump.



Lock and Unlock Keypad

The **UNO+** device can lock and unlock the keypad by pressing and holding the **Silence Alert Button** for 5 seconds at any time during therapy.



Sleep/Power Saving Mode

The device goes into **Sleep/Power-Saving** mode after 2 minutes of inactivity. Press and hold any button for 1 second, after which the device will no longer be in **Power-Saving** mode. Press the **desired button** to change the therapy settings.



Pressure/Mode Setting Selection

Select the desired pressure setting by pressing the Therapy button. If the LED light on the left side is **ON**, the device is in Continuous mode. If the LED on the right side is **ON**, the device is in Variable mode. To change from Continuous to Variable mode, press the **Pressure/Mode Button**.



Therapy Pressure Selection

Select therapy pressure of 40mmHg, 60mmHg, 80mmHg, 100mmHg or 125mmHg by pressing the **Pressure Settings Button**. The green LED light next to the button will indicate the selected pressure. Variable Mode options: 125mmHg/30mmHg, 100mmHg/30mmHg, 80mmHg/30mmHg. To select preprogrammed variable pressure, press the Therapy Mode button and select variable. Then press the pressure setting button and select 125mmHg, 100mmHg, or 80mmHg. The pressure will cycle every 3 minutes.



Troubleshooting Alerts

If the UNO+ device begins to alert, refer to the Troubleshooting Guide

THINGS TO KNOW

Battery Charge

- Only use the USB charger and USB cable provided with the UNO+ device.
- CAUTION: Do not place the device in its carrying case during battery charging.
- Connect the USB cable to the charger head. Make sure you have fully connected the USB cable.
- Connect the USB charger to a nearby outlet. Plug the USB cable into the lower part of the UNO+ NPWT device (slot below the Power Button).
- Let the device charge for approximately 5 hours. The LED will blink every 3 seconds, which indicates that the device is being charged. After 5 hours, the battery will be fully charged, and the LED next to the Low Battery Icon will remain steadily ON as long as the charger is plugged in. Once the charger is disconnected from the device, the LED will turn OFF

Insert Canister

Attach the canister to the back of the device. Slide in the canister and make sure it clicks to confirm that the canister is secure and tightly held.

Canister Removal and Replacement

Turn the **UNO+** device **OFF**. Slide the dressing tubing clamp close to where the tubing connects to the canister. Close the clamp. Disconnect the port tubing from the canister. Remove from the case if in use. Depress the tab on the two sides of the canister to remove the used canister. Install a new canister. Reattach the dressing tubing to the canister tubing port. Return the device to the carrying case, if in use. Release the tubing clamp. Turn the UNO+ device **ON**.

QUESTIONS?

Call: +1.800.208.2025 or +1.516.487.8787

Email: info@genadyne.com | www.genadyne.com

the power to heal