

Highly cost effective customer options

Whether you chose to purchase or rent, XLR8 provides THE most cost effective choice available, guaranteed!

- Ultra low cost Pump rental or purchase
- Similar priced consumables

Don't be scalped on price from more expensive competitors!

XLR8+

A more powerful pump which easily handles the largest of wounds for extended periods of time. The added advantage of being WiFi enabled which provides:

- GPS asset management
- Real time monitoring of patient's usage
- Improved user compliance
- Real time alarm notification and troubleshooting during an event
- Photo capturing and uploading to the patient's file on the webserver
- Remote control for caregiver to change settings in real time over the internet
- Software updates
- Ability to program hours of usage
- Cost analysis information
- Data analysis



XLR8 and XLR8+



the power to **heal**
in the palm of your hand

Genadyne Biotechnologies
16 Midland Avenue
Hicksville
New York
11801
USA

Phone: +1 516 487 8787
Toll free: +1 800 208 2025
Fax: +1 877 487 7878
Email: info@genadyne.com
DMR-06-094 Rev A

XLR8 Components



Green Foam

- A fully reticulated polyurethane foam made with polyether resin to allow better hydrolytic stability, providing increased performance whilst under pressure
- Duo-Decahedron structure gives the cells three-dimensional skeletal strands, providing enhanced filtering qualities
- The combination of cell structure and composition gives the foam a non-adherent property, for easier removal during dressing changes and almost eliminating pain
- Promotes perfusion and assists in tissue granulation
- Fraying is significantly reduced, preventing stray pieces from being left in the wound bed



PVA White Foam

- Polyvinyl Alcohol Foam moistened with sterile water
- High Tensile strength
- Comes in Small, Medium and Large size kits



Silicone Port Pad

- Softest patient interface – enhanced patient comfort less chance of pressure points
- Film based secure system – for easy application and dressing security
- Integrated clips – quick and easy changes and facilitates patient showering



Bridge Dressing

Spiro Foam

- Size: 104 +/- 1.5 cm uncoiled length x 1.9 +/- .3 cm width x 1.6 +/- 0.15 cm height
- Comes in Small and Medium size kits



XLR8 Drape

- 3 step application guide – easy application
- Pre-marked cutting guides – for easy and accurate application
- Clear visibility – for easy alignment over the foam
- Enhanced film security – ability to stay in place for extended periods of time

Pump

- Lightweight (585grams) – easy to maneuver and carry
- Whisper quiet operation – patients get better quality sleep
- 12 hour lithium battery with 3 hour recharge – easy to manage, minimal interruptions to care
- Supports multiple languages – clinical staff, patients and caregivers can understand the functions of the pump
- 5 advanced alarms – leak, blockage, canister full, low battery and critical battery
- Easy and intuitive programming and maintenance
- Clear and visible display – Pressure, system and battery status
- Simple lockable controls – guarantees no patient/ caregiver identifiable interference with care
- Comes with carrying case and strap plus bed holder



Canister

- 5 capacity canisters 200ml, 400ml, 600ml, 800ml and 1100ml – allows all exudate levels to be managed
- All canisters fit onto the pump and fit into the carry case/bed holder
- Slip and clip attachment to pump – easy to change
- Mix and match capability – use any canister with any dressing kit – any wound size/exudate level can be managed optimally

Y connector

Allows multiple wound sites to be treated with one pump – enhanced cost effectiveness



Dressing kits

- Available in various formats:
- Thick Green Foam – Small (S), Medium (M), Large (L) and Extra Large (XL)
 - Gauze (AMD) Kits – S, M and L
 - Silver Foam Kits – S, M and L
 - PVA (White Foam) – S, M and L
 - Abdominal Kit
 - Spiro Foam Kit
 - Bridge Dressing Kit

XLR8 NPWT

The XLR8's advanced variable pressure capabilities allow different negative pressures to be used without the downsides of Intermittent Therapy.

By maintaining negative pressure throughout the cycle, it avoids losing an airtight seal and displacement of the dressing, which would allow exudate to pool in the wound bed.


This gradual “wave” is gentler to the patient, increasing blood flow to the peri-wound, and preventing ischemia. The changes in pressure are imperceptible to the patient, increasing patient comfort and compliance.

Without the sudden peaks in pressure from Intermittent Therapy, in-growth of the granulation tissue to the foam dressing is virtually eliminated which makes dressing changes less stressful and less painful to the patient.

Advanced Intermittent Variable Negative Pressure Wound Therapy

Intermittent pressure levels of other brands = 125-0mmHg (spike)

Variable pressure levels of XLR8 = 80-30mmHg (gradual)

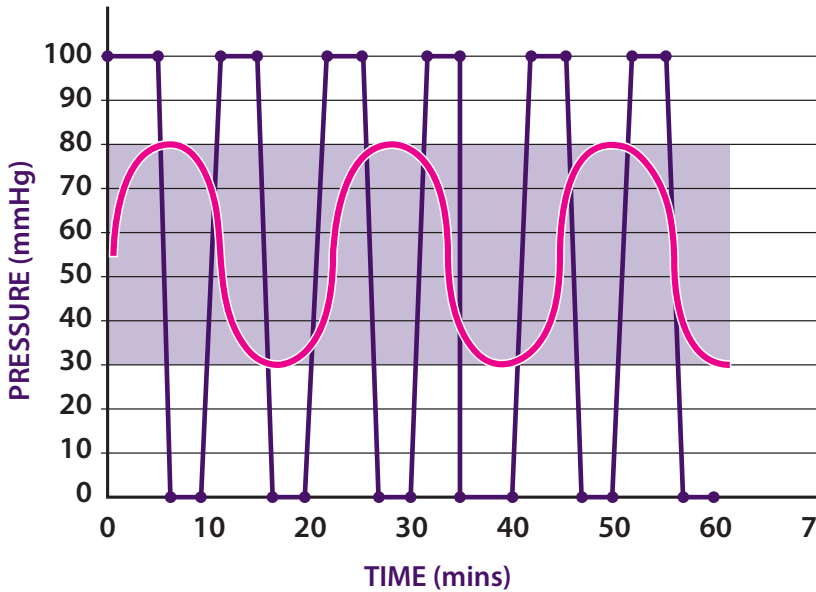

Benefits of XLR8

- Completely customizable high and low pressure settings
- Always maintains negative pressure at wound site
- Eliminates discomfort caused by pressure spikes
- Based on clinical needs, pressure can be set above or below capillary pressure of 32mmHg